

Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India P: 9210 903 903, +91 129 4272727, © 93111 24302 E: info@prestogroup.com



Faridabad
Sonipat
Kolkata
Mumbai
Pune
Ahmedabad
Chennai
Bangalore
Hyderabad

www.prestogroup.com





Model No. PAC-0123

Autoclave is a consistent tool used to sterilize various medical tools in microbiology labs, pharmaceutical, veterinary science, dentistry, plastics, production lineups etc. PRESTO Autoclave is made of high grade stainless steel material with radial locking system.

- Microprocessor based display for accurate test results
- Highly accurate test results under pressure saturated steam
- Radial Locking system
- High Grade stainless steel construction
- Autocut off Safety Device for temperature & pressure overshoot
- Bright LCD display
- · Feather touch controls
- Pressure regulating valve for maintaining pressure inside test chamber

Autoclave is used to create sterilization requiring elevated temperature and pressure different from the ambient air pressure. Additional Options like hysteresis, purging, proportional band available for optimization process.

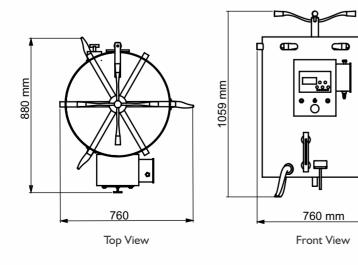

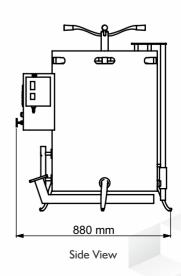



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India P: 9210 903 903, +91 129 4272727, S 93111 24302 E: info@prestogroup.com



• Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad • Chennai • Bangalore • Hyderabad

www.prestogroup.com

Key Specifications:

Display	LCD (Digital)
Pressure Gauge Range	0-30 psi, 2.1 Kg/cm²
Timer	Upto 999.9 min
Power	220V , Single phase, 50 Hz
Temperature Range	working range 121°C / Max. 125°C
Test Bucket Size	Internal Dia 12" x Height 17"
Water Requirement	5-8 Ltr.

Highlights:

Settings	Chamber temperature, Sterilization time, Dry time, High Temperature
Auto cutoff Facility	Yes
Material	High Grade Stainless Steel (304 Grade)
Dimensions	760x880x1060 mm
Temperature Control	PID
Water Level Indicator	Yes

Optional Features:

Available in various capacity as per requirement.

Thank you customers for choosing us as your partners in growth!































OUR OTHER PRODUCTS:-













Top Load Tester

Wall Thickness Gauge

Gloss Meter Triangle

Eccentricity Tester

Moisture Analyzer MB 45

Preform Neck Cutter

^{*}All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.