



M 95

Moisture Analyzer



Combining Speed and Precision in One Moisture Analyzer – Test More Samples in Less Time


The M95's halogen heating system distributes heat uniformly over the sample and speeds up the drying process — delivering fast and precise measurements. Designed for tool-free cleaning, the M95's components can be easily removed for quick and easy maintenance. Operating the M95 is intuitive with the icon-driven menu navigation on the touchscreen display.

Standard Features:

- **Increase Efficiency in your Lab with Fast Heating Halogen Technology**
The MB95's precisely controlled halogen heating dries samples quickly and ensures uniform heating to yield fast, precise and repeatable results of 0.01% readability. Four drying profiles and five switch-off criteria help to perform customized sample tests.
- **Define your Drying Parameters and Store Up to 20 Methods for Easy Recall**
MB95 can help you store up to 20 methods and 200 total results, for statistical analysis of measurements.
- **Effortless Maintenance with the Tool-Free Cleaning Design**
Designed with effortless maintenance in mind, the MB95 features a removable pan handler and heating chamber components that make cleaning quick and easy—no tools needed.
- **Touchscreen Icon-Based Menu Navigation Guides Users through Every Step**
Icon-based menu navigation on the easy-to-read touchscreen display makes operating the M95 intuitive. Perform routine sample testing in three easy steps —no training required.



Moisture Analyzer

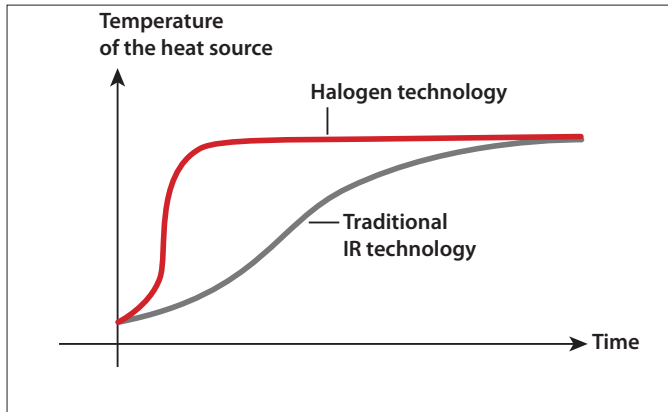

Exceptional Precision

The M95 provides precise results for more confident testing.

- **Fast Halogen Heating Technology**
Halogen heating dries samples in seconds and performs up to 40% faster than traditional infrared methods.
- **Uniform Drying for Precise, Repeatable Results**
The reflective sample chamber ensures uniform heating of each sample to yield reliable results of 0.01% readability.
- **Flexible and Powerful**
Features a full operating temperature range of 40– 200 °C in one degree increments.
- **Four Drying Profiles to Ensure Accuracy**
Customize heating of specific samples with temperature control for one-step drying—set temperature and time and begin!
Fast: Quickest possible time to temperature
Standard: Minimal overshoot while achieving temperatures quickly
Ramp: Controlled ramping to temperatures
Step: Set three temperatures at three time intervals
- **Shut-Off Criteria for Automatic Test Completion**
Customize your own test endpoint criteria to automatically end the test once moisture loss is no longer detected. Choose from five criteria— Manual, Timed, A30, A60, A90.

Advanced Functionality

Test a wide variety of materials for a range of applications with M95

- **20 Methods Library Makes Testing Faster**
Store methods for up to 20 different samples; methods can be easily recalled and transferred to other units for quick setup.
- **Store Up to 200 Results for All Methods for Statistical Analysis**
The M95 will store 200 results which can be printed or exported via USB memory stick for further analysis
- **RS232 and USB Ports**
Quick connections to a printer, PC or flash drive





Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

- Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
- Chennai • Bangalore • Hyderabad

www.prestogroup.com



M 95

Easy Tool-Free Cleaning

The M95 features a patented pan handler and removable parts for convenient maintenance.

- **Stainless Steel Drying Pan**
Wipe down quickly
- **Removable Heating Chamber Components**
No tools required—simply remove glass and tray to clean



Easy to Use and Maintain

The M95 features a wide, bright touchscreen display and intuitive menu navigation for fast, accurate moisture analysis.

- **Graphic Display with Touchscreen**
Navigate the device easily with an intuitive user interface on a bright, color, and easy to read graphic display
- **Icon-Driven Menu Navigation**
Easy to follow — no training required
- **Four Mechanical Keys for Fast One-Step Operation**
Tare, Start/Stop, Print, On/Off



www.prestogroup.com



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India

P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com

• Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
• Chennai • Bangalore • Hyderabad

www.prestogroup.com



M 95 Moisture Analyzer

Model	M 95
Capacity	90 g
Readability	0.01%/0.001 g
Repeatability (Std Dev) (g)	0.1% (3 g sample) 0.018% (10 g sample)
Moisture Range	0.01% to 100% (0.01% to 1000% for regain mode)
Heating Element	Halogen
Drying Programs	Standard, Fast, Ramp, Step
Temp Range	40° - 200° C
Switch-off Criteria	Manual, Timed, Auto (30, 60, 90 seconds)
Adjustment	External adjustment mass - 50 g
Power	100 - 120 VAC 5A 50/60 Hz 200 - 240 VAC 2.5A 50/60 Hz
Operating Temperature Range	5° to 40°C
Display Type	4.3" QVGA, TFT touch screen
Display Results	%moisture, %solids, %regain, time, temperature, weight, method name, drying curve and statistics
Pan Size (mm)	Ø 90
Interface	RS232, USB host, USB device
Adjustable Feet and Level	Yes
Dimensions (W x H x D) (mm)	210 x 180 x 350
Net Weight	5.23 kg
Shipping Weight	8.35 kg

Compliance

- **Product Safety:** IEC/EN 61010-1, -2-010; CAN/CSA C22.2 61010-1, -2-010; UL 61010-1, -2-010
- **Electromagnetic Compatibility:** IEC/EN 61326-1 Class B, Basic environments; FCC Part 15 Class A; ICES-003 Class A
- **Compliance Marks:** CE, CSA, RCM, WEEE

Other Standard Features & Equipment


ABS top housing, stainless steel pan support, stainless steel pan handler, in-use cover, Five Switch-off criteria (Manual, timed, A30, A60 or A90), Four heating profiles (standard, fast, ramp and step), storage up to 20 methods and 200 test results for each method, user guide, 40–200 °C heating temperature range (1 °C increments), 14 operational languages (display only, no key in)

Accessories

Temperature Calibration Kit, MB series11113857
In-Use-Cover.....30284478
Filter, Glass Fiber, Set (200)80850087
Cage, Sample80252477
Anti-Theft Device80850043

Pan Holder.....30246947
Pan, Reusable, Set (3), 7mm, MB series80252478
Pan, Reusable, Set (3), 14mm, MB series.....80252479
Aluminum pan packed with 50 pcs30585411
Scoop30284477

Outline Dimensions :

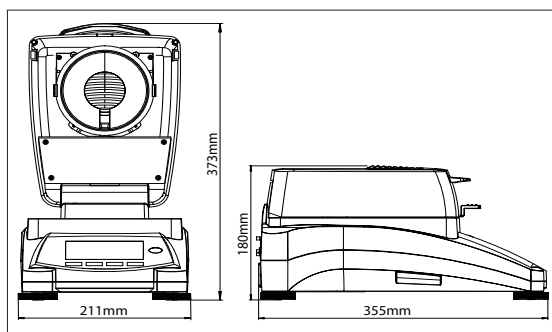

www.prestogroup.com