



Cyclic Salt Spray Test Chamber 2.0



Model No. CCT-450 2.0

Registered Design No. 244580

Presto Cyclic Corrosion Salt Spray Chamber - Worlds most Interactive and versatile corrosion test systems range available, synthesizing salt spray with other environmental conditions such as fixed humidity, air drying and wetting / condensation.


Presto Cyclic Corrosion Salt Spray Chamber models will perform the neutral salt spray tests, along with in addition these CCT models are loaded with extra features which enables users to perform intense 'modified'. In this the conventional salt spray is often programmed with multiple climate, in a multi-stage cycle, to intensify the test. Some Test Cycles Names - Salt spray and Wetting / Condensation Humidity (SWAAT) or Salt spray and Air drying (PROHESION).

Complaint to ISO 16701, ISO 11997-1, VDA 233-102, JASO M609, ASTM G85

- HMI based Touchscreen Model
- Relative Humidity can be monitored on HMI Screen.
- Cyclic Test mode enabled to set parameters according to test sample
- Corrosion Resistant PT-100 sensor
- Solid state digital Programmable Temperature Controller
- Program profile for selecting test parameters and performing the test.
- Easy Data Management. User can create PID settings, Air Saturator settings, test temperature, test and air purge time etc.
- Current test running status display
- Robust Construction with Fiber reinforced body(Triple wall)
- Ease of use of sample placements inside the test chamber
- Pneumatically operated canopy
- Air purge after test for fog settlement.
- Reservoir Tank inside sponge filter with low solution level alarm.
- 5 point Auto Tuning feature selection
- Real time vs temperature graph
- Data recording with temperature data reading
- USB option available (Pen drive is not a part of supply)
- User can control PID settings through Advanced Hysteresis setting
- Angle of Canopy maintained in order to restrict flow of droplets directly on the test sample.
- Low compressed air alarm (Smart device introduce)
- Low NACL solution alarm.
- Over / under voltage protection
- Individual specimen real time monitoring with identification.
- Maintain fault history (Easy to Diagnost the exact issue)

New Innovative Add On features in Presto Cyclic Test Chamber that make it Class Apart.

- Ethernet Connectivity- to connect with other devices using LAN and internet.
- Also to connect any computer or other electronic devices to its network as long as the device has an Ethernet adapter or network card.
- Built-in VNC server that allows remote monitoring and operating of the HMI functions via VNC Client APP (Windows, iOS, Android)
- Web Monitoring through which allows direct monitoring of register data via web -page, and requires no additional software I installation
- Alarm Notification-Automatically sends out alarm notification emails to logged-in recipients when alarms occur and supports the Secure Sockets Layer (SSL) protocol to ensure safe data transmission.
- Option to interchange HMI fonts between HINDI and ENGLISH for Operator convenience to operate and access database using Hindi Language.

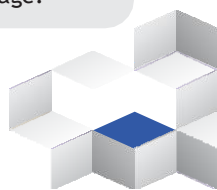



Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
 P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com



- Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
- Chennai • Bangalore • Hyderabad

www.prestogroup.com

Key Specifications:

Display	HMI Based touchscreen 7" TFT LCD Color screen
Humidity (RH Range)	30-95 % RH \pm 5% RH
Dry Cycle Temperature Range	25 to 60 °C \pm 3°C
Salt Cycle Temperature Range	35 °C \pm 2°C
Chamber Temperature Least Count	0.1°C
Air Saturator Temperature Range	40°C - 55°C
Air Saturator Temperature Least Count	0.1°C
Air Saturator Temperature Repeatability	\pm 2°C
Power	220, Single phase, 50 Hz, 15 A
pH value of Salt Solution	6.5 to 7.2 pH/3-3.3 pH
Temperature Control	Inbuilt PID Temperature Controller
Timer(Hour meter)	6 digits
Test Air Pressure	0.7Kg/cm ² - 1.2Kg/cm ² *
Air Regulator	0 to 30 psi
NACL Solution Tank Capacity	60 Ltr.

Highlights:

Air Saturator water Level low alarm	Yes
Fog Collection Funnel	1
Air Purge	Auto/Manual
Air Saturator Auto Refilling	Yes
Material	Fiber Glass
Over Temperature Alarm	Yes
Reservoir float valve	Yes

* Compressed Air requirement : Approx 7 cfm (Air Compressor not part of supply)
 Main air line pressure regulator with filter in customer scope.

Optional:

- Online pH meter.
- Table top pH meter.

MODEL	Internal Size (in mm)
Cyclic 450 Liters	1000 x 700 x 750
Cyclic 1000 liters	1600 x 700 x 800



Control Panel with Ethernet port

Check for options



Pneumatic Canopy lift system

* available with only 250 liter model



Automatic Leveler For Air Saturator



Solution Level Indicator with Alarm



Specimen Positioning System




Smart Monitoring device under & over voltage protection



Air Monitoring device with Low Air Alarm



Inbuilt Humidity Source

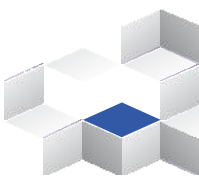





Presto Stantest Pvt. Ltd.

I-42, DLF Industrial Area Phase-1, Delhi Mathura Road, Faridabad 121003, Haryana, India
P : 9210 903 903, +91 129 4272727, 93111 24302 E : info@prestogroup.com



- Faridabad • Sonipat • Kolkata • Mumbai • Pune • Ahmedabad
- Chennai • Bangalore • Hyderabad


www.prestogroup.com

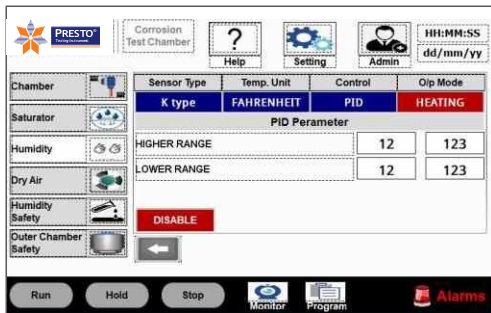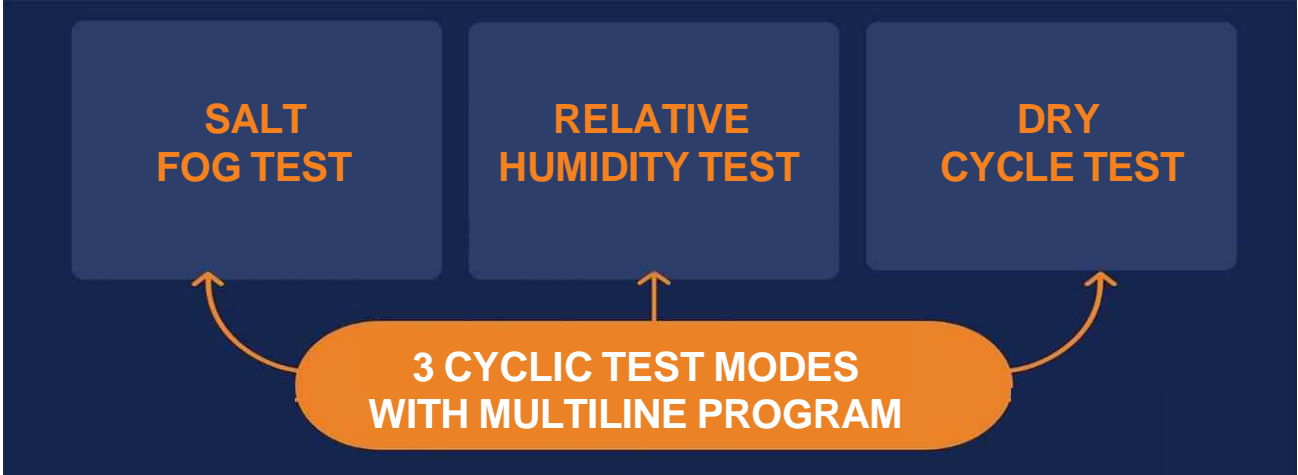


ADMIN SETTING



INDIVIDUAL TEST

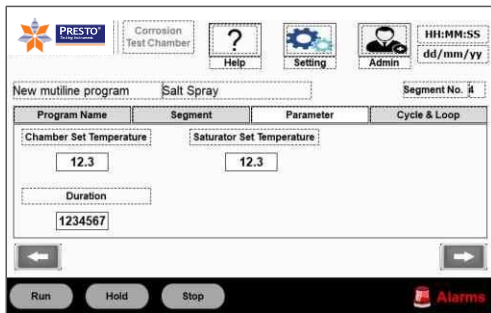

MULTI LINE PROGRAM

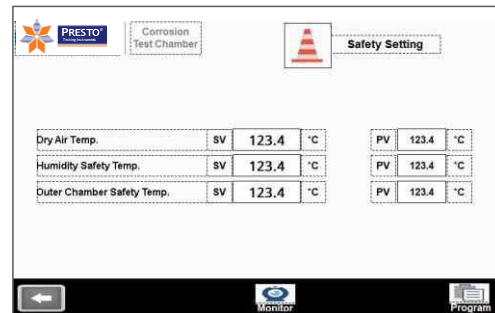

SAFETY SETTING

Thank you customers for choosing us as your partners in growth !



*All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.

