


Scuff Resistance Tester



Model No. PSR-292

Scuff Resistance Tester is a unique test instrument to determine the color transfer from printed or coated material against rubbing. It checks the color transfer when one surface is rubbed against other when subjected to continuous rubbing under specified load. PRESTO Scuff Resistance Tester is designed as per standards ASTM F2497-05, BS 3110:1959

- Motorised Weight arrangement on sample for noise free working
- Highly accurate test results under uniform load and specified movement.
- Uniform load by weights on test sample for evenness in testing
- Provision of Set number of counts through Digital Preset Counter
- Firm hold of test specimen by circular clips/clamps
- Motorised Geared system for scuff movement
- Digital Pre-settable counter upto 9999 counts
- Highly accurate and precise results under uniform load and rubbing clamp movement
- Accurate rpm of upper rotating discs using synchronous motor
- User friendly electronic controls and ergonomically designed feature
- Adjustable footrest for leveling.

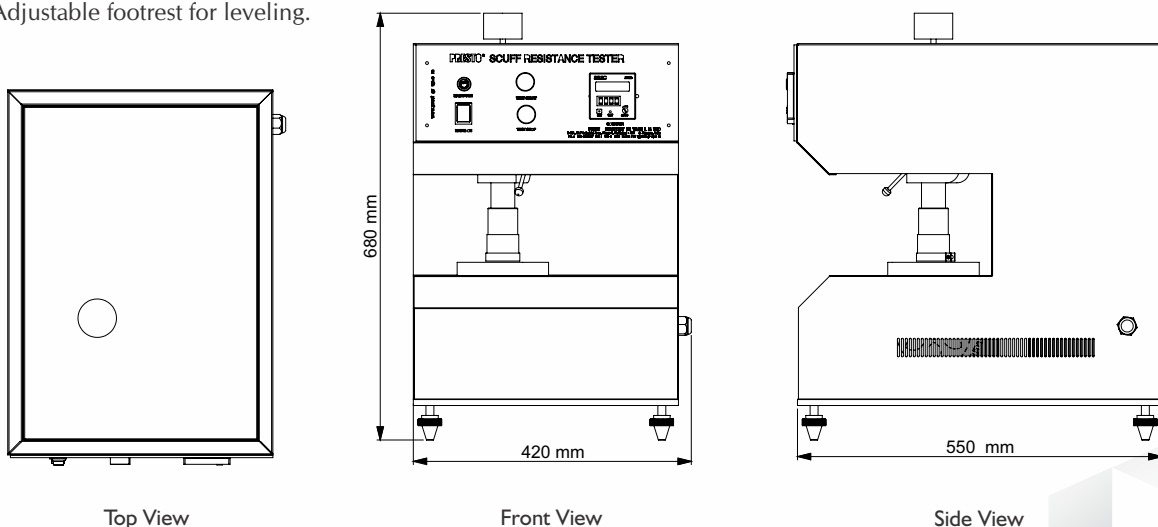

It is designed and produced according to ASTM & BS Standard and measures the number of double strokes need to cause Fuzz (appearance of loss of fuzz) and Peel (actual rupture of the surface).

Key Specifications:

Test Weight	2 psi
Counter	4- digit Digital Preset type
Max.range of counter	9999
Least Count of counter	1
Rotational Speed	60 rpm \pm 2 rpm
Lower Motor capacity	¼ HP
Upper Motor capacity	10 kg/cm ² torque motor
Diameter of Upper Clamp	48 \pm 2 mm
Diameter of Lower Clamp	120 \pm 2 mm
Power	220V, 50Hz, single phase

Optional:

- Dust removing blower system on request.
- Weights above 2 psi available on request.

Highlights:

Display	Digital (LED)
Safety Lock	Yes (with upper clamp)
Sample Preparation templates	3
Material	Mild Steel
Finish	Powder coated Havel Gray & Blue combination finish and bright chrome / zinc plating for corrosion resistant finish
Dimensions	420x550x680 mm

Thank you customers for choosing us as your partners in growth !



*All the above logos are the sole property of their respective owners and are used purely for depiction purposes only.

OUR OTHER PRODUCTS :-



Head Office : Presto Stantest Private Limited

I-42, DLF Industrial Area Phase-I, Delhi Mathura Road, Faridabad 121003, Haryana, India
 P : +91 129 4272727 E : info@prestogroup.com

Branches :

Mumbai

503-504B, Fair link Centre, Opp. Fame Adlab
 Behind Mongini's Off Link Road
 Andheri (West), Mumbai 400058
 P : +91 022 42957052, 42957053
 E : west@prestogroup.com

Kolkata

93, 3rd Floor, Karnani Estate
 209, A.J.C. Bose Road
 Kolkata-700017. West Bengal
 P : +91 033 40612012
 E : east@prestogroup.com

Chennai

A-3, First Floor, Sri Sathya Sai Building
 755, Poonamallee High Road, Kilpauk
 Chennai- 600 010. Tamil Nadu
 P : +91 044 43865260, 43838372
 E : south@prestogroup.com

